


HORIZON™ MANUAL RESUSCITATOR

LATEX FREE | SINGLE PATIENT USE | DISPOSABLE

The Horizon™ Manual Resuscitator allows you to deliver a high standard of patient care while effectively managing healthcare costs. The resuscitators come standard with a bacterial/viral filter, oxygen reservoir and mask. Optional safety components and comfort features allow you to create a custom resuscitator to meet the needs of each patient and of your facility's budget.

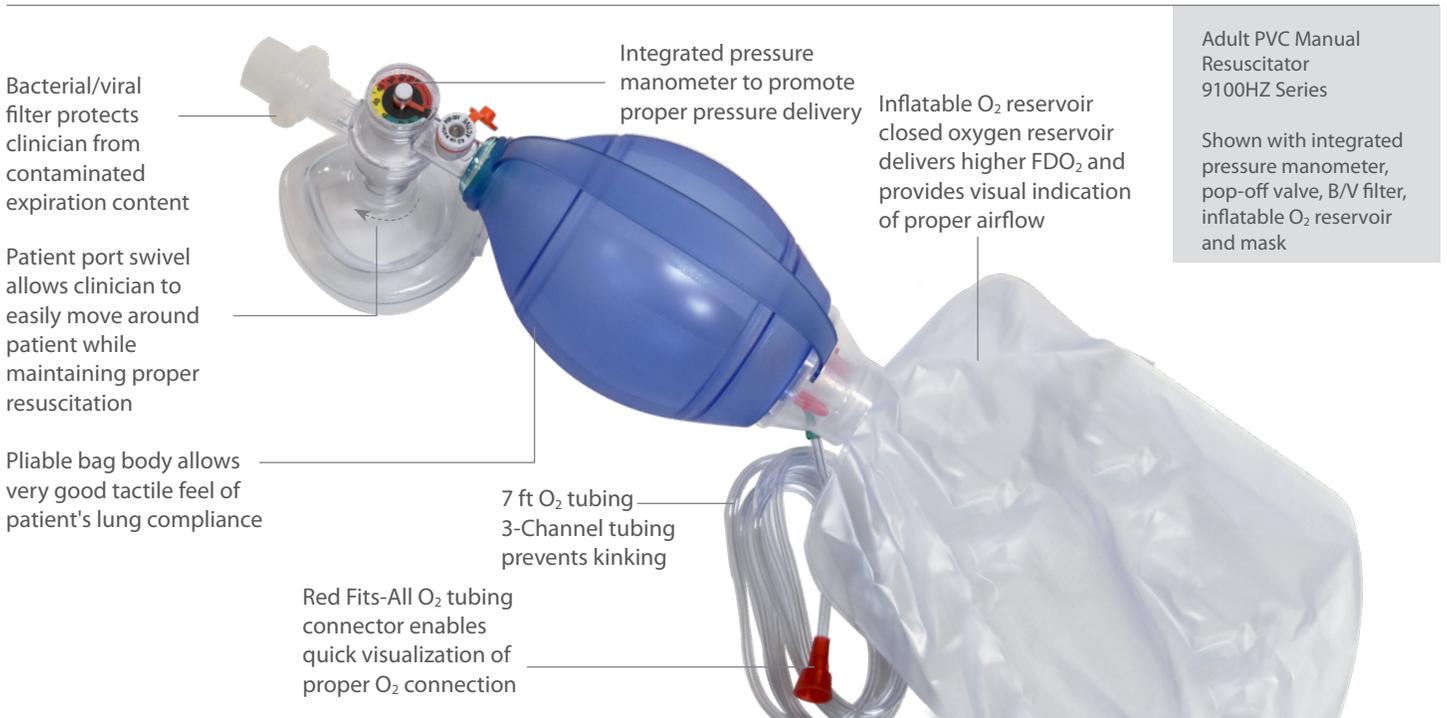

HORIZON™ MANUAL RESUSCITATOR

- Single patient use resuscitation/ventilation bag available in four sizes/volumes
- Available with or without integrated pressure manometer
- Integrated color-coded pressure manometer effectively promotes proper pressure delivery (< 20 cm H₂O); the tri-color fields alert to the risk level of delivered pressure
- Accurate pressure delivery helps reduce gastric insufflation, thus preventing vomiting, aspiration and resulting pneumonia
- Face mask proportionately sized to resuscitation bag size
- Select from a variety of components to assemble a custom resuscitation bag to meet your unique patient needs (see opposite side)

PVC MANUAL RESUSCITATOR

ITEM	SIZE	VOLUME	STROKE VOLUME	PK
6100HZ	Infant	300 mL	140-170 mL	10
8100HZ	Pediatric	750 mL	≥ 450 mL	10
9100HZ	Adult	1625 mL	≥ 600 mL	10

SILICONE MANUAL RESUSCITATOR

ITEM	SIZE	VOLUME	STROKE VOLUME	PK
6000HZ	Infant	300 mL	140-170 mL	10
7000HZ	Child	550 mL	340-360 mL	10
8000HZ	Pediatric	1000 mL	590-610 mL	10
9000HZ	Adult	1500 mL	690-710 mL	10

HORIZON™ MANUAL RESUSCITATOR

CUSTOMIZE YOUR RESUSCITATION / VENTILATION BAG

Choose from the following components to create a custom resuscitation/ventilation bag to best meet your needs:

OXYGEN RESERVOIR



INFLATABLE O₂ BAG

Closed system oxygen reservoir delivers higher FDO₂ and provides visual indication of proper airflow



CORRUGATED TUBING

For use with secondary oxygen supply; available in 24 inch and 48 inch lengths



COLLAPSIBLE/EXTENDABLE O₂ ACCUMULATOR


Available only for 9000HZ and 9100HZ resuscitators

CONVENIENCE COMPONENTS



HAND STRAP

Maintain a positive grip with a hand strap
**Only available on inflatable bag reservoir*

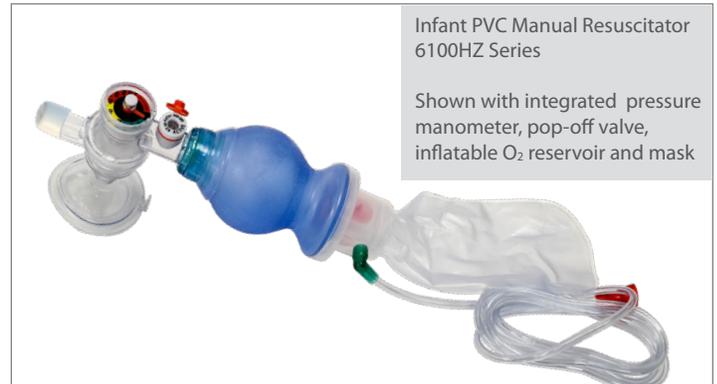

Infant PVC Manual Resuscitator
6100HZ Series

Shown with integrated pressure manometer, pop-off valve, inflatable O₂ reservoir and mask

SAFETY COMPONENTS



INTEGRATED PRESSURE MANOMETER

Ventlab's innovative color-coded dial clearly denotes pressure delivery:

- Green** - Target pressure level <20 cm H₂O
- Yellow** - Potentially dangerous pressure level
- Red** - Dangerous pressure level



POP-OFF VALVES

40 cm H₂O & 25 cm H₂O
Allows clinician to use an appropriate pop-off level



FILTERS

A bacterial/viral filter provides effective protection for the healthcare provider during expiration



PEEP VALVES

Adjustable PEEP valves;
5-20 cm H₂O