


COMBIGUARD II™ + COMBIPORT®

LATEX FREE | STERILE | SINGLE PATIENT USE | DISPOSABLE

The Combiguard II™ and Combiport® are easy to use, highly effective wound irrigation tools. The Combiguard II Irrigation Splash Guard covers the wound and helps prevent splash-back. The Combiport attaches to a standard saline bag to deliver a high volume of fluid rapidly. Available sterile and non-sterile.

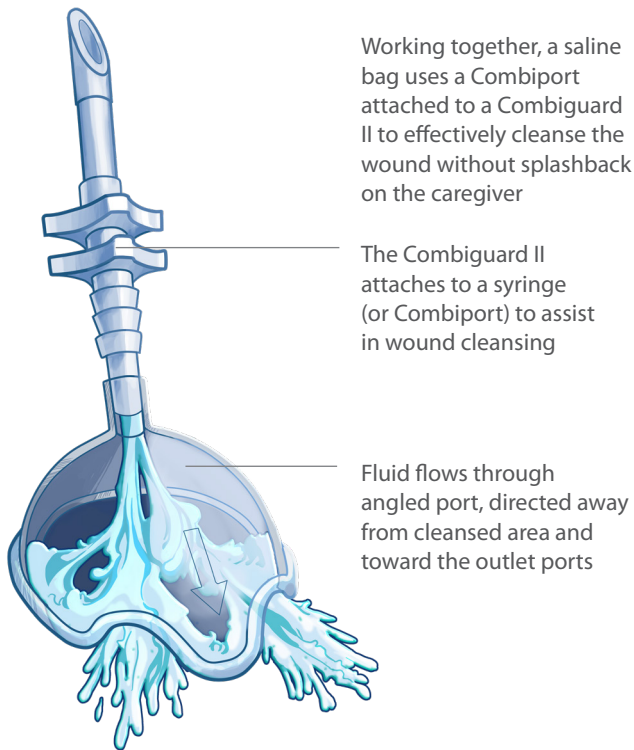

COMBIPORT® WOUND IRRIGATION DEVICE

- Attaches to a standard saline bag to administer a high volume of fluid in a short period of time
- Irrigation helps remove foreign bodies, reduce pathogen numbers and moisten wound tissue
- Patented spike attaches to standard 500 mL or 1000 mL saline bag
- For use only with Combiguard Irrigation Splash Guard (not sold separately)
- Removable cap for easy stop/start during cleansing
- Easy to assemble and use

STERILE

ITEM	DESCRIPTION	PK
10222	Combiguard II Irrigation Splash Guard	50
10100	Combiguard II + Combiport	50

NON-STERILE

ITEM	DESCRIPTION	PK
10222NS	Combiguard II Irrigation Splash Guard Non-Sterile	500

COMBIGUARD II™ IRRIGATION SPLASH GUARD

- Egg-shaped cup covers wound site to prevent irrigation splash back, helping to protect clinician
- Dual outlet ports enable contaminated fluids to safely escape
- Fluid delivery inlet delivers a stream equivalent to an 18 g needle
- Angled delivery inlet directs fluid away from cleansed area toward outlet ports
- Rounded edges provide a smooth surface against the skin
- Fits on 30-40 cc luer-tipped syringe
- Can be combined with Combiport Wound Irrigation Device