

# ARTERIAL BLOOD GAS (ABG) SAMPLING

LATEX FREE | SINGLE USE | DISPOSABLE


Arterial Blood Gas (ABG) sampling is an invasive procedure that can result in accidental needlestick injury and exposure to bloodborne pathogens. The Pulset™ ABG Syringe features SafetyTip™ Needlestick Protection to help prevent accidental needlestick injuries in the care setting. The SafetyTip™ is engineered as an integral part of the needle to comply with all FDA, OSHA and NIOSH guidelines. The Pulset™ has other additional features making it one of the easiest and quickest ABG syringes to use in healthcare today.



Puncture Kit



A-Line Kit

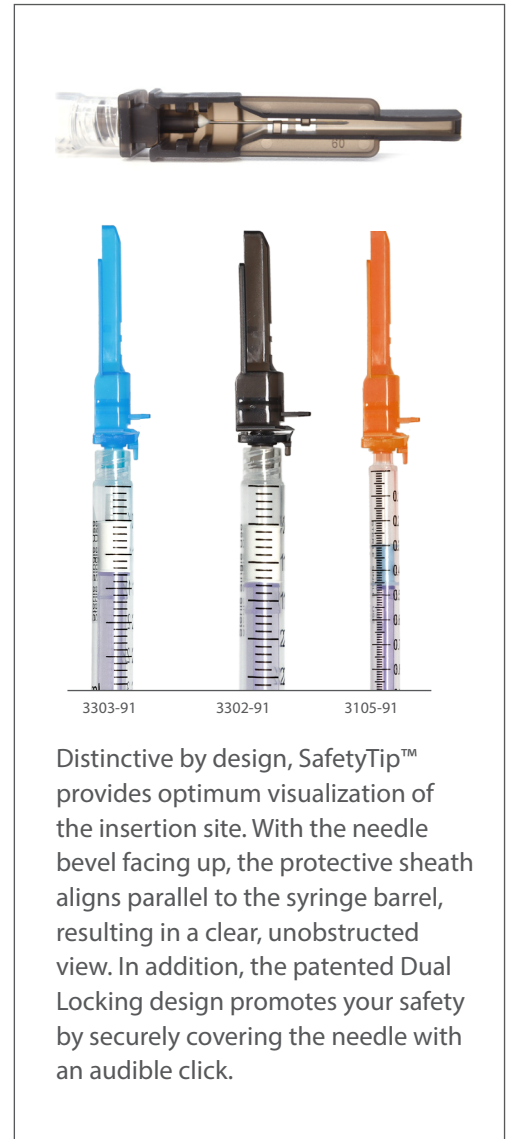

Distinctive by design, SafetyTip™ provides optimum visualization of the insertion site. With the needle bevel facing up, the protective sheath aligns parallel to the syringe barrel, resulting in a clear, unobstructed view. In addition, the patented Dual Locking design promotes your safety by securely covering the needle with an audible click.

## PULSET™ ADVANTAGES

### SHARPS PROTECTION DEVICES (SAFETY TIP™)

- The SafetyTip comes attached to the needle with bevel in the ready/up position
- The SafetyTip device folds back along the syringe for optimal visualization of the insertion site
- Single hand activation

### VENTED PLUNGER

- Eliminates the need for manual clearance of air bubbles
- Allows for fewer steps between needle insertion and sending the sample to the lab

### DRY BALANCED LITHIUM-ZINC HEPARIN

- Eliminates calcium electrolyte errors
- Easily mixes with blood to prevent clotting

### THIN-WALLED NEEDLE TECHNOLOGY

- Flows as quickly as standard needles
- Smaller puncture site means less pain for the patient

1cc Pulset Syringe, Vented

ITEM	DESCRIPTION	ALCOHOL PREP	IODOPHAR PREP	GAUZE	ICE BAG	BANDAGE	Q-CORK	PATIENT LABEL	QTY
3103-91*	23g x 1" Needle with SafetyTip, 25u BH								100
3105-91*	25g x 5/8" Needle with SafetyTip, 25u BH								100
3165-91	Modified Kit, 25g x 5/8" Needle with SafetyTip, 25u BH	X	X	X	X	X		X	100
3173-91	Full Kit: 25g x 5/8" Needle 23g x 1" Sterile Unattached Needle, 25u BH	X	X	X	X	X	X	X	100

3cc Pulset Syringe, Vented

ITEM	DESCRIPTION	ALCOHOL PREP	IODOPHAR PREP	GAUZE	ICE BAG	BANDAGE	Q-CORK	PATIENT LABEL	QTY
3302-91*	22g x 1" Needle with Safety Tip, 25u BH								100
3303-91*	23g x 1" Needle with Safety Tip, 25u BH								100
3362-91	Modified Kit, 22g x 1" Needle with SafetyTip, 25u BH	X	X	X	X	X		X	100
3362-91S	Modified Kit, 22g x 1" Needle Luer Slip with Safety Tip, 25u BH	X	X	X	X	X		X	100
3363-91	Modified Kit, 23g x 1" Needle with SafetyTip, 25u BH	X	X	X	X	X		X	100
3373-91	Full Kit: 23g x 1" Needle with SafetyTip, 22g x 1.5" Sterile Unattached Needle with Safety Tip, 25u BH	X	X	X	X	X		X	100

Cord Gas Kit

ITEM	DESCRIPTION	ALCOHOL PREP	IODOPHAR PREP	GAUZE	ICE BAG	BANDAGE	Q-CORK	PATIENT LABEL	QTY
3182-91	Cord Gas Kit, (2) 22g x 1" Syringes with Needles with Safety Tip, 25u BH				X			X (2)	100

With Blunt Tip Cannula

ITEM	DESCRIPTION	QTY
3100-30*	A-Line 1cc Syringe Blunt Tip Cannula 25u BH	100
3300-30*	A-Line 3cc Syringe Blunt Tip, 25u BH	100

A-Line

ITEM	DESCRIPTION	QTY
3100-20*	A-Line 1cc Syringe 3cc Waste Syringe 25u BH	100
3100-25*	A-Line 1cc Syringe 25u BH	100
3300-21*	A-Line 3cc Syringe 5cc-6cc Waste Syringe 25u BH	100
3300-25*	A-Line 3cc Syringe 25u BH	100

\*Items include CE Mark.